



A GUIDE TO FOLDING WIPES

1

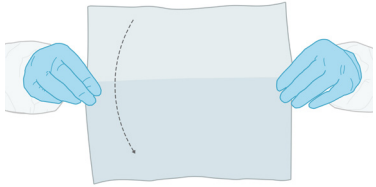


Hold the wipe in mid-air. Take the wipe by its edges and fold it in half. The wipe should not be placed on any surface during this folding action

2

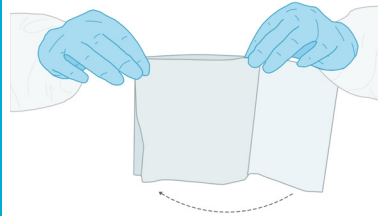

Fold the wipe in half again, creating four cleaning faces - two outside and two inside

3

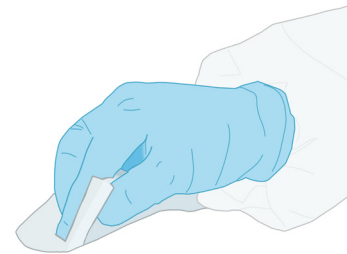


Wipe in one direction with slightly overlapping strokes from the cleanest areas to the dirtiest areas, lifting the wipe at the end of each stroke

4




Each face of the wipe will be used for one straight stroke. Turn the wipe over and make a second stroke

5

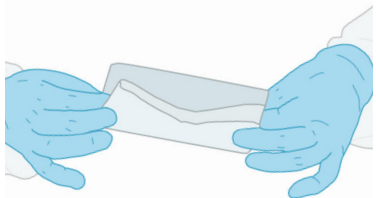


Refold the wipe as shown in step 2 to enclose the two contaminated faces and expose the two clean faces. Repeat step 3 & 4 using the two clean faces

6

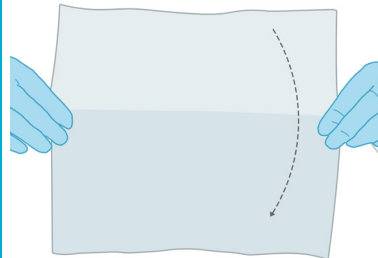


The wipe may be completely unfolded, reversed, and refolded beginning at step 1 to enclose the contaminated faces and expose the four remaining clean faces

7

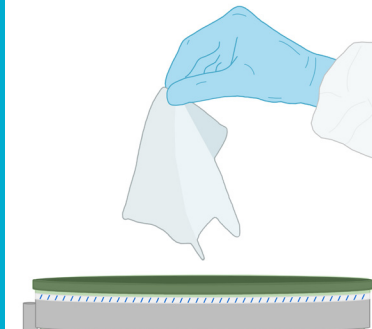

Once all 8 sides have been contaminated, dispose of the wipe according to site protocol

8

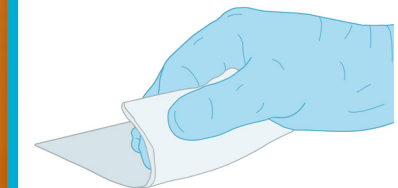

Continue with a new clean wipe, following all prior steps until the full surface area is clean